



Camino[®] Intracranial Pressure Monitor

Empowering clinicians to achieve early and accurate detection of unpredictable events

The Camino is an advanced, touch-screen monitoring system that captures real-time data at the time of catheter insertion and trend data over a five day period. Camino is a compact and portable device that provides tools for continuously determining and monitoring ICP and ICT intracranially.




Advanced and versatile features:

- Tunneled or bolted catheter options
- Intuitive, easy-to-use menu system puts five days of patient information at your fingertips with rapid data export to USB
- Widescreen display ensures that critical data can be read from 10 feet away and up to a 30 degree offset from the monitor
- Ergonomic handle and convenient size for bedside pole attachment

Versatility in draining and monitoring with a family of bolted and tunneled catheters

Natus Camino catheters measure their respective values at the tip of the catheter, providing continual ICP measurements. This design eliminates the need for a fluid filled system to communicate pressure to an external transducer.

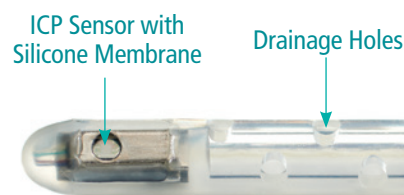

Camino Flex Catheter Tip

Camino Family of Catheters

Camino Flex - Therapeutic ventricular drainage of CSF while performing continuous ICP monitoring with the Camino Monitor

B-Series - Parenchymal bolt catheters for measuring ICP and ICT

HM-Series - Ventricular bolt catheters for measuring ICP and ICT

1104G - Subdural, post craniotomy catheter for measuring ICP

1104L - Parenchymal bolt catheter that secures through the Licox® bolted system for measuring ICP

Detailed information about the Natus Camino Catheters can be found in our product catalog.

Key Features of Camino catheters:

- Easy and reliable way to measure Intracranial Pressure and Brain Temperature
- Camino bolt is quick and easy to place – Parenchymal and ventricular options available
- Accurate monitoring at the source
- Continuous parenchymal or ventricular monitoring
- Ease of use - catheters can be easily connected and disconnected from the monitor
- Proven technology



Camino Flex Catheter

Service

Excellence in customer and technical service is a core value of Natus Neuro.

Essential Support Elements:

- Accessible and effective Technical Support
- Excellent replacement unit or spare part availability
- Extended warranty and service coverage programs
- Comprehensive, flexible customer training courses

To learn more about Natus Neuro products contact your local distributor or sales representative.

US Customers Call: **1-800-356-0007**

International Customers Call: **+1-608-829-8500**

natus
neuro

Natus Neuro

3150 Pleasant View Road
Middleton, WI 53562 USA

Tel: 1-800-356-0007

1-608-829-8500

Fax: 1-608-829-8709

www.natus.com